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ENGLISH

ELASTIC THUMB SPLINT WITH SPIRAL STAY

REF 03-1-004

DESCRIPTION

Elastic thumb splint featuring a flexible spiral stay for thumb support. Equipped with adjustable hook-and-loop (Velcro) fastening straps. Helps restrict thumb movement while providing support and stabilization. Washable. Fits either the left or right hand. Color: Beige.

INDICATIONS

- Conservative treatment of thumb tendon and ligament injuries.
- Thumb joint arthritis.
- Thumb dislocations and ligament injuries.
- Suitable for first metacarpal fractures.
- De Quervain's syndrome.

WARNINGS

- Do not use this product without prior advice from your healthcare professional or qualified therapist.
- Do not use the product if the size is incorrect. Selecting the correct size is essential for proper fit, optimal performance, and user safety.
- Always apply the product to clean and dry skin.
- Do not apply the product over areas with skin lesions, bleeding disorders, swelling, or varicose veins without consulting your physician or healthcare professional.
- Discontinue use if discomfort, skin irritation, or pain occurs.
- Do not use the product if it is deformed or damaged.
- Do not attempt to repair or modify the product yourself.

SAFETY INSTRUCTIONS

- Read these instructions carefully before using the product.
- Always follow the instructions for use to ensure safe and effective application.
- Inspect the product before each use for signs of wear, tears, or deformation. Do not use the product if any damage is detected.
- Incorrect use or failure to follow these instructions may result in complications or injury.
- This product is intended for use by a single patient only and should be used according to the instructions of a physician or physiotherapist.
- If pain occurs, persists, or worsens, discontinue use immediately and consult your physician or physiotherapist.
- The manufacturer is liable only for manufacturing defects, provided that the product has been used in accordance with these instructions.
- The manufacturer's liability and war-


ranty apply only to products purchased through authorized distributors or specialized retailers.

CLEANING & MAINTENANCE

- Clean the product regularly, at least once a month or more frequently if necessary, to maintain it in optimal condition.
- Hand wash with lukewarm soapy water at a temperature of up to 30°C.
- Do not iron, bleach, dry clean, or tumble dry.
- Store the product in a clean, dry place and avoid exposure to extreme temperatures.

INSTRUCTIONS FOR USE

1. Apply the splint by positioning the section with the flexible spiral stay along the thumb.
2. Secure the thumb strap using the hook-and-loop (Velcro) closure.
3. Secure the wrist strap using the hook-and-loop (Velcro) closure.
4. Check for proper fit.

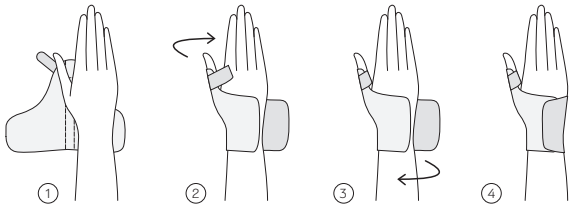ΜΕΓΕΘΗ / SIZES	SMALL	MEDIUM	LARGE	X-LARGE
Περίμετρος καρπού Wrist circumference	16-18 cm	18-20 cm	20-22 cm	22-24 cm

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